

# TruPort

Flexibility has never been easier



Enhancing outcomes for  
patients and their caregivers:

**Hill-Rom**

# Ready for every assignment

## TruPort: functional solutions for greater flexibility and improved ergonomics

Medical workstations in ICUs and operating rooms are subject to new challenges and modifications almost daily. An increasing number of devices and instruments are improving medical possibilities but require extremely flexible supply systems: not just adjustable storage areas or extendable connections for gas, electricity and data, but improved safety and ergonomics too. It is imperative that everything is in the right place.

The highly flexible, ceiling-mounted supply system, TruPort, by Trumpf Medical satisfies all end user requirements. Its modular design allows you to integrate new components whenever necessary or indeed to completely reconfigure the supply units. This versatility makes TruPort one of the most efficient and future-proof supply systems in the world.



## Your benefits at a glance

- **Future-proof**

The flexible architecture of TruPort makes it easy to switch innovative components and accessories including the supply modules for gas, electricity and data

- **Adaptability**

The TruPort support heads can be easily and quickly equipped with supply modules, different components and accessories on all four sides and along their entire length

- **Ergonomics**

The color-coded brake design, intelligent work-station components and optional motorized lifting mechanisms provide for a comfortable and intuitive system

- **Intelligence**

The optional MPC Rail is the central point of control on the support head and is also used to attach innovative accessories and supply them with power

- **Performance**

Large coverage area, high load capacity and almost unlimited configuration options make a broad range of applications possible


- **Design**

Closed, flush-fitting surfaces help achieve hygiene requirements and make the system aesthetically pleasing



# A masterpiece of modularity

As a revolutionary, modular system, TruPort meets all user requirements. In addition to the supply modules, the support head can be equipped with different shelves and drawers to suit individual requirements. At any time.

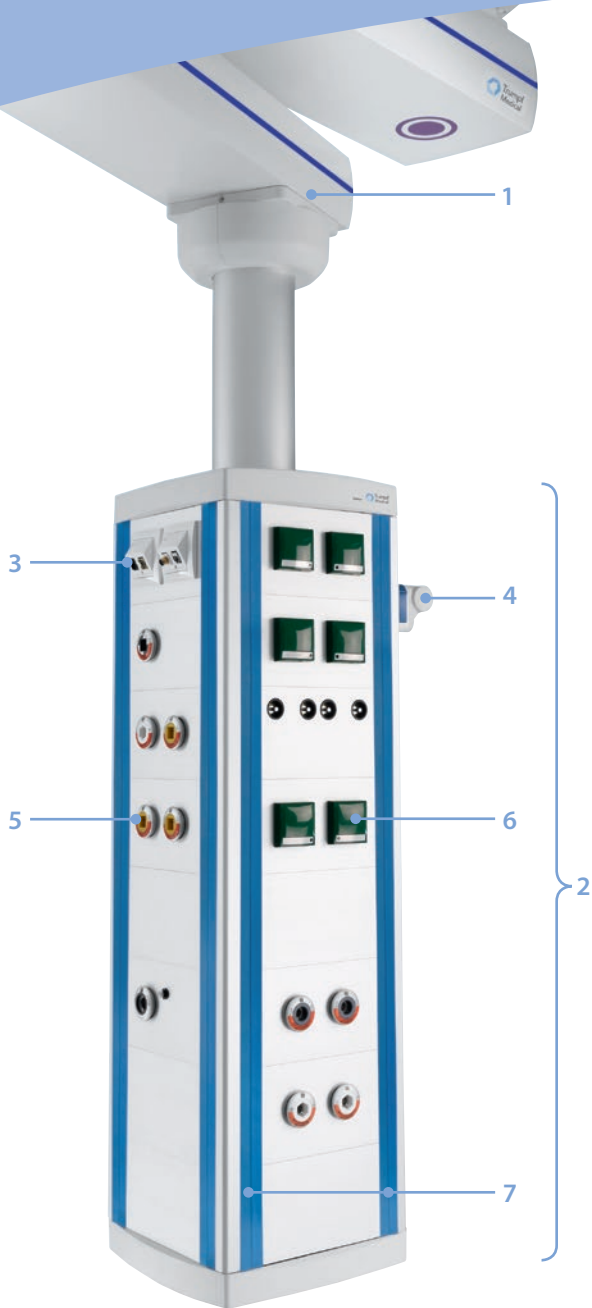

- 1 Support arm, dual-arm system
- 2 Support head size L; length 39"/1,000 mm
- 3 Supply module for data
- 4 LED lights (only MPC Rail)
- 5 Supply module for gas
- 6 Supply module for electricity
- 7 Rail

## Configure with complete freedom



Medium-load shelf, maximum of two drawers, with optional control handle



Heavy-load shelf, maximum of four drawers, with optional control handle



The heavy-load shelf with keyboard drawer makes it easy to document in the room; its handle-free design also makes it easy to clean

# Achieving new standards in usability



Conversion of a supply module for electricity

## A freely configurable central supply unit

Every TruPort support head can be equipped with different supply modules required, from gas to electricity and data. Trained personnel can make additions or changes easily without the need for common tools, even after initial installation.



A click mechanism provides for straightforward, tool-free fixing of workstation components such as shelves on every side and at all heights of the support head

## The right components for every workstation

The click mechanism allows different components to be fixed to the integrated M Rail at all heights and on every side, without the need for tools. A wide range of components are available for individual workstation design. Shelves of different sizes for medium and heavy loads provide space for every type of medical-technical device. Different drawer designs support ongoing processes and make it easier to work.

# Flexible for greater ease of use

## Flexible and intelligent: TruPort with MPC Rail

Special electronic workstation components can be attached to the optional MPC Rail which uses low voltage power. Function commands are forwarded via a BUS system, providing for new possibilities for workstations that are more innovative and ergonomic than ever before.




Click mechanism for comfortable, tool-free insertion of pivoting LED lamp



Ease of use thanks to remote control of MPC Rail – for raising the anaesthesia device and operating the brake



MPC Rail: intelligent with 24 V rail and integrated BUS system

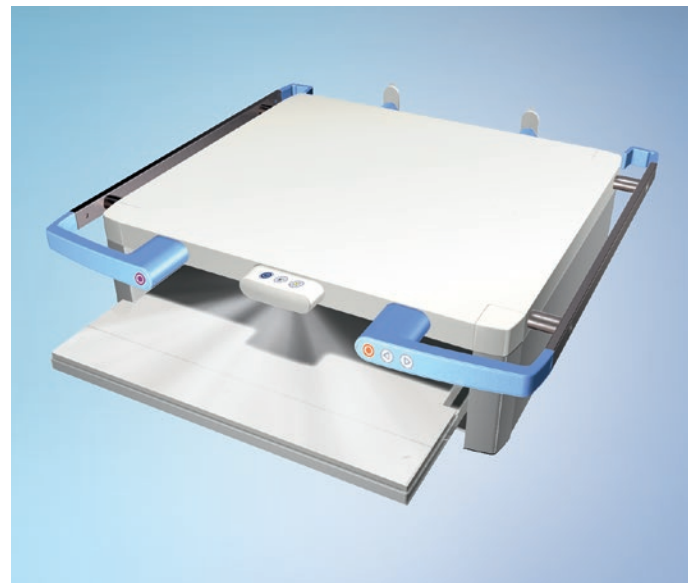

Of special benefit in the ICU: drawer with AmbientLine desk illumination

## The elements of the TruPort modular system and their features at a glance

	standard	add. options with MPC Rail
<b>Support arms</b>		
Support arm reach	<b>single arm system:</b> up to 51"/1,300 mm <b>dual arm system:</b> up to 126"/3,200 mm	
Brake system	friction pneumatic	electro-pneumatic brake system electro-magnetic brake system <sup>[1]</sup>
Ceiling attachment	suspended ceiling clearance up to 98"/2,500 mm seismic approved	
Height adjustment	●	
Payload	up to 750 lbs/340 kg	
<b>Support head</b>		
Dimensions	M (12" x 9"/310 x 220 mm) L (12" x 12"/310 x 310 mm)	integration of motorized lifting system
Lengths	up to 69"/1,750 mm	
<b>Supply modules</b>		
Data module reorganization	●	
Gas module reorganization	●	
Electric power module reorganization	●	
<b>Workplace components</b>		
Toolless mounting with click mechanism	●	electronic components
AmbientLine modules	AmbientLine Floor, AmbientLine Ceiling	AmbientLine Desk
<b>Shelves and drawers</b>		
Load capacity	<b>medium load shelf:</b> 110 lbs/50 kg <b>heavy load shelf:</b> 176 lbs/80 kg	
Operating element on shelves	for medium load shelf	for heavy load shelf
Shelf dimensions (WxD)	<b>medium load shelf:</b> 21", 25" x 19"/ 530, 630 x 480 mm <b>heavy load shelf:</b> 17" <sup>[2]</sup> , 21", 25", 29" x 19"/ 430 <sup>[2]</sup> , 530, 630, 730 x 480 mm <b>depth extension kit</b> 3"/70 mm available	
Number of drawers per shelf	Max. 4	

● = available

<sup>[1]</sup> Country-specific <sup>[2]</sup> 430 mm (17") width also available in 260 mm (10") depth

Subject to change without notice. The information provided in our quotation and order confirmation is decisive.



For more about TruPort and its technical features and configuration options, go to [www.trumpfmedical.com/en/truport](http://www.trumpfmedical.com/en/truport)

Hill-Rom is a leading global medical technology company with more than 10,000 employees in over 100 countries. We partner with health care providers by focusing on patient care solutions that improve clinical and economic outcomes in five core areas: **Advancing Mobility, Wound Care and Prevention, Clinical Workflow, Surgical Safety and Efficiency, and Respiratory Health**. Hill-Rom people, programs, and product brands work towards one mission: **Every day, around the world, we enhance outcomes for our patients and their caregivers.**

Trumpf Medical, part of Hill-Rom, is distinguished by high-quality German engineering standards and offers innovative products to improve efficiency and safety in the OR, ICU, and in other clinical environments throughout the care sector. With our customers' requirements as our benchmark and innovation as the foundation of our success, Trumpf Medical delivers total solutions to fit your clinical care needs.

This document is destined solely for use by healthcare professionals. Medical devices shown in this brochure are intended for use with patients in departments of healthcare establishments.

These products are regulated health care products which, where required by applicable regulations, bear a CE mark. Hill-Rom recommends you carefully read the detailed instructions for safe and proper use included in the documentation accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2015 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.  
Doc. No: 2027337, 30 September 2015

Not all products/options are available in all countries.  
For further information about our products or services, please contact your local Trumpf Medical representative or visit our webpage:

[www.trumpfmedical.com](http://www.trumpfmedical.com)

Enhancing outcomes for  
patients and their caregivers:

**Hill-Rom**